

SURGICAL

“MOLDED” N95 RESPIRATOR

NIOSH Surgical Rated Approved
 TC-84A-9290 (Regular), TC-84A-9292 (Small)

Part Number: PT-N95C/CS-02S
Description: Molded N95 rated respirator with nylon head bands
Location: Manufactured in Ohio, USA
Brand: Patriot Mask



KEY FEATURES:

- NIOSH approved N95 rating
- Meets NIOSH inhalation and exhalation resistance testing procedures: CFR 84.172 .
- ≥95% BFE (Bacterial Filtration Efficiency) according to ASTM F2101.
- ≥95% Filter Efficiency according to ICS Laboratories Inc.
- Pre-formed for secure fit.
- Nose foam
- Moldable nose clip will help provide custom fit and secure seal.
- Stapled headbands
- Available in Small and Regular
- Latex free: Yes
- Model: Molded Cup

DO NOT USE:

- For atmospheres less than 19.5% oxygen
- For gases and Vapors, including oil aerosols, etc.
- For repeated use following NIOSH guidelines
- If modified or used for its unintended purpose
- Expiration date of 3 years after manufacture date

MATERIAL COMPOSITION:

- Outer Layer: Spun bound cloth/non-woven fabric.
- Straps - Polyisoprene
- Staples - Ultrasonically welded
- Nose Clip - Aluminum
- Filter - Polypropylene
- Shell - Polyester
- Inner Layer: Needle punch cotton
- Nose Clip: 1 mm Thick Metal
- Nose Foam: Spandex filament sponge bar
- Headband: Nylon string
- Breathing valve: N/A

DO USE FOR:

- Medical clinics and operating rooms
- Laboratory procedures
- Infection control requirements
- Emergency planning and preparedness

SURGICAL

“MOLDED” N95 RESPIRATOR

NIOSH Surgical Rated Approved
 TC-84A-9290 (Regular), TC-84A-9292 (Small)


Part Number: PT-N95C/CS-02S
Description: Molded N95 rated respirator with nylon head bands
Location: Manufactured in Ohio, USA
Brand: Patriot Mask

PACKAGE INFORMATION

Description	Description	TYPE	PCS/Box	Box/Carton	Pcs/Pkg	Carton/Pallet	Total Quan
PT-N95C-02S	Regular	Surgical	25	20	500	30	15,000
PT-N95CS-02S	Small	Surgical	25	20	500	30	15,000

DIMENSIONS

REGULAR

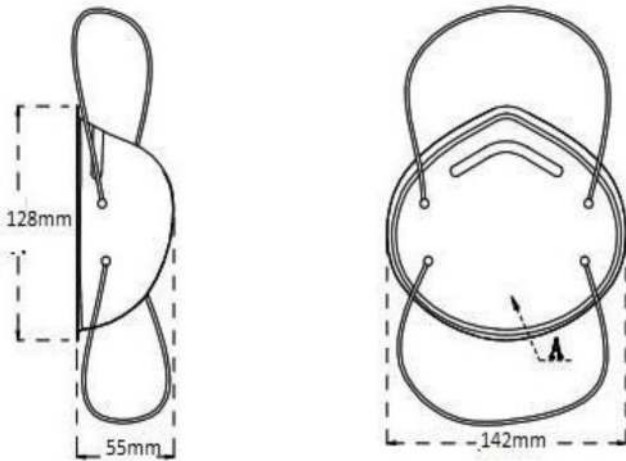

SMALL

